

Extract from the full product range

SAUER Leg bags

Inlet with push-on adapter, soft tube, 45 cm, adaptable length
The bags are available in boxes of 10 or 30 pieces.

Swing tap

Adapter	Capacity	Quantity	Order no.
smooth push-on adapter	0,6 l	10 pcs	710.1345-10
		30 pcs	710.1345-30
stepped push-on adapter	0,6 l	10 pcs	710.1445-10
		30 pcs	710.1445-30
smooth push-on adapter	1,0 l	10 pcs	711.1345-10
		30 pcs	711.1345-30

Sliding tap

smooth push-on adapter	0,6 l	10 pcs	710.2345-10
		30 pcs	710.2345-30
stepped push-on adapter	0,6 l	10 pcs	710.2445-10*
		30 pcs	710.2445-30*
smooth push-on adapter	1,0 l	10 pcs	711.2345-10
		30 pcs	711.2345-30

Twist tap

stepped push-on adapter	1,0 l	10 pcs	711.3445-10*
		30 pcs	711.3445-30*

* available also as sterile version

SAUER Sport Bags

with or without fixed adapter, capacity 0,5 litres
Die Beutel sind in Packungen mit 10 oder 30 Stück erhältlich.

Swing tap

Adapter	Inlet	Quantity	Order no.
fixed adapter, smooth	diagonal inserted tube	10 pcs	721.1100-10
		30 pcs	721.1100-30

Sliding tap

fixed adapter, smooth	diagonal inserted tube	10 pcs	721.2100-10
		30 pcs	721.2100-30

SAUER Rectangular Bag

With fixed adapter, capacity 0,6 litres
The bags are available in boxes of 10 or 30 pieces.

Sliding tap

Adapter	Inlet tube	Quantity	Order no.
fixed adapter, stepped	32 cm, smooth, nicht kürzbar	10 St.	732.2232-10
		30 St.	732.2232-30



MANFRED SAUER GMBH

Hersteller von Medizinprodukten

SAUER Leg Bags

... for the able-bodied and for wheelchair users



Different shapes of bags with a wide range of features



MADE IN GERMANY We develop and manufacture our products at Lobbach, Germany



MANFRED SAUER GMBH
Hersteller von Medizinprodukten

Neurott 7-15 • 74931 Lobbach • Germany
Tel. +49 6226 960 104 • Fax +49 6226 960 011
export@manfred-sauer.de • www.manfred-sauer.com

For further information about our wide range of urine bags please give us a call or send us an e-mail.

Tel. +49 6226 960 104
export@manfred-sauer.de

Urine collection bags

The urine collection bags have been developed and are produced by SAUER. They are manufactured on high-frequency welding machines at our premises at Lobbach. Their different shapes, sizes and various features make it possible to meet the particular needs and requirements of most users.



SAUER Leg Bags for able-bodied persons

(also for people in a wheelchair)

- Rectangular shape; straight inlet, capacity 0.6 litres or 1.0 litres; optionally with backing of viscose flock on one side
- Inlet: push-on adapter, smooth or stepped; soft tube, 45 cm, can be shortened; non-return valve
- Outlet: swing tap, twist tap or sliding tap




SAUER Sport Bags

- Straight or diagonal inlet; capacity 0.5 litres
- Inlet: without adapter or with fixed adapter, smooth or stepped; non-return valve
- Outlet: swing tap, twist tap or sliding tap

SAUER Rectangular Bag

- Capacity 0.6 litres
- Inlet: fixed adapter, stepped, oblique inserted tube, fixed length, 8 cm or 32 cm, non-return valve
- Outlet: sliding tap

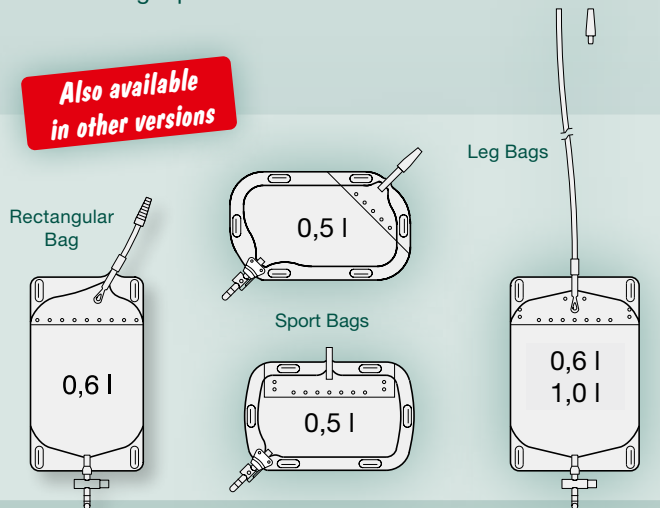

Properties

- The Connector**
with or without adapter, allows a simple and safe connection for all catheters and all urinary sheaths:
 - smooth adapter: connects to urinary sheaths
 - stepped adapter: connects to catheters
- The Inlet Tube**
soft version and in different lengths – can also be shortened to adapt to the body height.
- The Welded Valve**
allows for a smooth flow of urine, even with leg bags which do not lie flat.
- The Non-Return Valve**
prevents urine backflow e.g. during physiotherapy sessions.
- The Eyelets**
allow the straps to be fastened according to individual need, thus giving a more secure fit.
- Der Auslaufhahn**
that allows draining the bag in a clean fashion:
 - Swing tap: Easy to handle, also with impeded dexterity of the hands
 - Twist tap: Compact and discreet, particularly for small bags, connects to a night bag easily
 - Sliding tap: Can be handled with only one hand; the bag can be worn on the right or on the left leg
- The Welding**
The welded valve, the non-return valve and the outlet tap are welded securely to the film. The double welding enables the bag to bear high strain.
- The Backing of Viscose Flock**
To avoid direct contact between the film of the bag and the skin, all urine collection bags (except the night bags) can be furnished with a backing of viscose flock on one side.

100% testing

To increase security for the user, every bag is *individually* tested for leaks and functional capability

