

Extract from the full product range

Urine Leg Bags for Wheelchair-users

Capacity 1.3 Liter

The bags are available in boxes of 10 or 30 pieces.

Swing tap

	Adapter	Inlet tube	Qty	Order no.
1.3 litres	without adapter	20 cm hard adaptable length	10 pcs	701.1020-10
			30 pcs	701.1020-30
	smooth adapter fixed	12 cm soft	10 pcs	701.1112-10*
			30 pcs	701.1112-30*

Sliding tap

1.3 litres	stepped adapter fixed	22 cm soft	10 pcs	701.2222-10*
			30 pcs	701.2222-30*
	stepped adapter push-on	35 cm soft adaptable length	10 pcs	701.2435-10*
			30 pcs	701.2435-30*

Twist tap

1.3 litres	stepped adapter push-on	35 cm soft adaptable length	10 pcs	701.3435-10
			30 pcs	701.3435-30

* available also as sterile version

Capacity 0.7 litres

The bags are available in boxes of 10 or 30 pieces.

Swing tap

	Adapter	Inlet tube	Qty	Order no.
0.7 litres	without adapter	20 cm hard adaptable length	10 pcs	702.1020-10
			30 pcs	702.1020-30
	stepped adapter fixed	12 cm soft	10 pcs	702.1212-10
			30 pcs	702.1212-30

Twist tap

0.7 litres	stepped adapter push-on	35 cm soft adaptable length	10 pcs	702.3435-10
			30 pcs	702.3435-30




MANFRED SAUER GMBH

Hersteller von Medizinprodukten

SAUER Bendi Bag

... particularly suited for wheelchair users

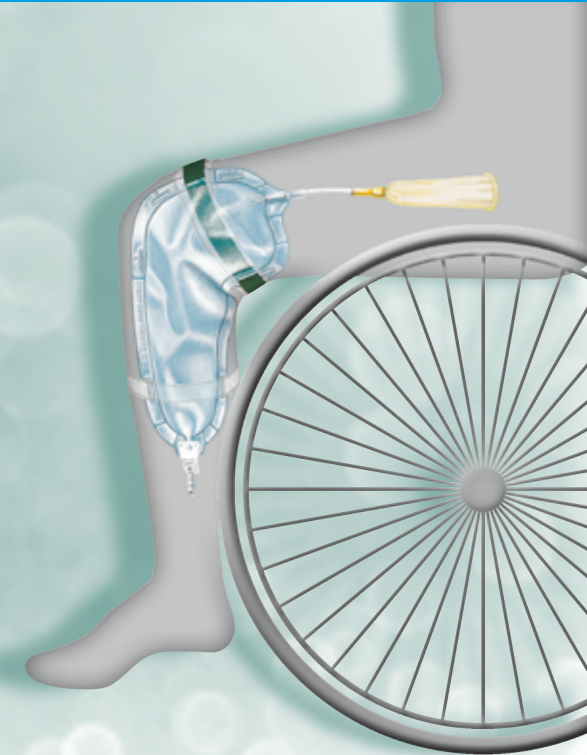

Urine collection bags,
available in different versions



MADE IN GERMANY We develop and manufacture our products at Lobbach, Germany



MANFRED SAUER GMBH
Hersteller von Medizinprodukten
Neurott 7 – 15 • 74931 Lobbach • Germany
Tel. +49 6226 960 104 • Fax +49 6226 960 011
export@manfred-sauer.de • www.manfred-sauer.com

For further information about our wide range of urine bags please give us a call or send us an e-mail.

Tel. +49 6226 960 104
export@manfred-sauer.de

Quality and wearing comfort

These bags with a leg hugging shape for fixation to the leg have been designed specifically for wheelchair users. They are adapted to the sitting position and offer a high degree of wearing comfort, flexibility and security. Due to their great variety of features, they fit every user according to his individual need.

SAUER Bendi Bag

High capacity

Available in 1.3l, the Bendi bag provides more independence on long journeys and at work without the worry of having to find accessible toilets. The smaller bags with a capacity of 0.7 l are better suited for children and youths.

Leg hugging shape

The angled inlet tube and the special shape are essential for a reliable use of these bags in the sitting position in a wheelchair, for comfortable and inconspicuous wear on the leg, a reliable urine outlet and easy emptying.

The textured film

This is very strong and has a slightly grained surface to avoid transpiration and prevent the film from sticking to the skin.


Also for able-bodied persons



Also available in other versions

100% testing

To increase security for the user, every bag is *individually* tested for leaks and functional capability.

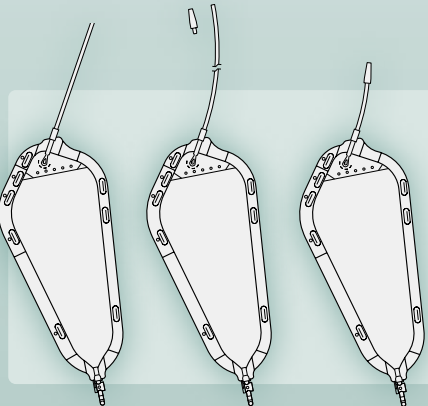

Properties

1 The Connector

With or without an adapter, allows a simple and safe connection for all catheters and all urinary sheaths:

- without adapter: connects directly to SAUER urinary sheaths
- smooth adapter: connects to urinary sheaths
- stepped adapter: connects to catheters



smooth adapter



stepped adapter



no adapter

2 The Inlet Tube

hard or soft version and in different lengths – can also be shortened to adapt to the body height.

3 The Welded Valve

allows for a smooth flow of urine, even with leg bags which do not lie flat.

4 The Non-return Valve

prevents urine backflow e.g. during physiotherapy sessions.

5 The Eyelets

allow the straps to be fastened according to individual need, thus giving a more secure fit.

6 The Outlet Tap

that allows draining the bag in a clean fashion:

- Swing tap: Easy to handle, also with impeded dexterity of the hands
- Twist tap: Compact and discreet, particularly for small bags, collects to a night bag easily.
- Sliding tap: Can be handled with only one hand; the bag can be worn on the right or on the left leg.

7 The Welding

The welded valve, the non-return valve and the outlet tap are welded securely to the film. The double welding enables the bag to bear high strain.



swing tap



swing tap reversed for left dressers



twist tap



sliding tap

